## SISSEL<sup>®</sup> HEATING PILLOW

---

---

**Instruction for use**  
**Gebrauchsanweisung**  
**Manuel d'utilisation**  
**Gebruiksaanwijzing**

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# SISSEL® Heating pillow

Thanks for choosing the SISSEL® heating pillow. To enjoy the comfort experience, please take a few minutes before starting operation of this pillow and read the following operation instruction carefully.

## IMPORTANT SAFETY INSTRUCTIONS RETAIN FOR FUTURE USE

The SISSEL® heating pillow is a high-quality product designed to last many years if it is properly used and cared for. Read the following instructions carefully to prevent damage or danger from incorrect use. Keep the instructions so that you can refer back to individual points when necessary.


If you let somebody else use the SISSEL® heating pillow, you should make sure the user has read the instruction before use.

### Items included in the package:

- 1 Heating pad
- 1 Inner cushion
- 1 Controller
- 1 Instruction for use

### Description:

1. Mains plug
2. Mains cable
3. Controller
4. Illuminated display for temperature settings
5. Sliding switch for ON/OFF and temperature settings
6. Plug-in coupling
7. Cover with integrated heating pad
8. Inner pad
9. Zip fastening

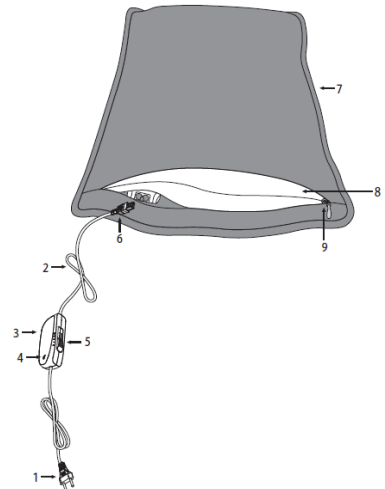**This appliance is not intended for medical use.**

### Warning:

Non-observance of the following information may result in personal injury (electric shock, skin burns, fire) or material damage. The following safety and hazard information is not only intended to protect your health and the health of others, it should also protect the product. For this reason, pay attention to these safety notes and include these instructions when handing over the product to others.

- Use the device only as an electric mattress heater.
- Only connect to the voltage shown on the rating plate.
- The appliance must not be used when folded or creased.
- Do not treat the heating pad as a sheet or truck it under the mattress.
- This appliance must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.

- Do not use the SISSEL® heating pillow on persons who cannot move without help, infants, persons who cannot feel heat, or animals.
- Persons who can feel heat to some extent may only use the SISSEL® heating pillow with supervision.
- Do not operate the SISSEL® heating pillow unsupervised.
- Caution! Never fall asleep while the SISSEL® heating pillow is in operation.
- Prolonged use can cause burns on the skin.
- Do not pierce the appliance with needles, pins or other sharp objects.
- Do not tug or twist the power lead.
- Do not expose the switch or power lead to moisture or liquids.
- In case of malfunction, do not use and contact the manufacturer.
- If the supply cord is damaged, it should be replaced by the distributor or its service agent or similar qualified persons in order to avoid a hazard.
- Do not use the appliance when it is damp or wet.
- The device may not be switched on using a time switch.
- The switch should not be placed on or beneath the appliance or be covered when the appliance is in use.
- Avoid sharp bending or bundling of the appliance.
- While the appliance is switched on:
  - Do not place any large or heavy objects on it such as suitcases, etc.,
  - Do not place any source of heat on it such as hot water bottles, heating pads or suchlike.
  - Do not allow the switch to be concealed under blankets, pillows or the like or lie on the SISSEL® heating pillow.
- Do not fold the SISSEL® heating pillow when it's hot.
- Examine the appliance frequently for signs of wear or damage. If there are such sign or if the appliance has been misused or does not work, do not use.
- The SISSEL® heating pillow is only to be used by children over the age of three under the supervision of parents or guardians. It should be ensured that they do not play with the appliance.
- Children below the age of 3 years are not allowed to use this device due to their inability to react to overheating.
- This appliance can be used by persons with reduced physical, sensory or mental capabilities if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the possible hazards involved.
- When the appliance is likely to be used for an extended period, the control should be adjusted to setting 1 for continuous use.
- When the appliance is not in use:
  - The plug should be removed from the mains and the appliance should be stored dry, without additional weight on it.
  - The connector should be disconnected from the appliance to avoid damage caused by twist or bend of the power cord.
- When not in use, store as follows:
  - When storing the appliance, allow it to cool down before folding.
  - Do not crease the appliance by placing items on top of it during storage.

**Switching on:**

Please always ensure that the SISSEL® heating pillow is not bundled or folded when it is in use.

The appliance may only be used with 230 Volts a.c.

Switch on the SISSEL® heating pillow using the hand-operated switch.

When switch on it has a 3 times beep –beep –beep sound to ensure the user has operated switch correctly.

Setting 0 = appliance off	Setting 2-3 = moderate warmth
Setting 1 = minimum warmth	Setting 4 = maximum warmth

**3 colors LED:**

- Yellow/Green light indicates warming up.
- Red light is a warning light: disconnect the SISSEL® heating pillow from power immediately!

**Quick heating:**

The SISSEL® heating pillow has a quick heating function that users can feel the warmth within 7-10 minutes when turn on the switch at highest setting.

**Timer function:**

Now you can experience the pleasant warmth. After about 90 minutes the SISSEL® heating pillow silently switches off automatically. In order to put the heating pillow into operation again you must first set the switch to position '0'. After about 5 seconds you can switch on the SISSEL® heating pillow again.

**Caution:**

The electronic components inside the switch may lead to a slight warming of the switch when the SISSEL® heating pillow is on. Therefore the switch must not be covered or placed on top of the SISSEL® heating pillow when it is in operation.

You should not rely on using your hands to gauge the SISSEL® heating pillow temperature because they are much less sensitive to heat than the surface of the rest of the human body. A temperature, which may feel tepid to your hand, may feel unpleasantly hot on your body.

**Cleaning and care:**

Please follow below symbols as indicated on the washing label.

## Heating Pad

1. Before cleaning the pad, first unplug it from the mains outlet, take out the inner cushion and let the appliance cool down.
2. The heating pad is equipped with a detachable cable. Disconnect the plug and remove the cable from the heating pad.
3. The switch or control unit must not be allowed to get wet during washing and during drying the cord must be positioned to ensure that water does not flow into the switch or control.
4. Please note that the heating pad should not be dry-cleaned, wrung out, machine-dried, pressed or ironed.
5. Set washing machine to a gentle washing cycle at 30° C. The heating pad should only be washed together with other textiles. Use a washing powder or liquid for delicate fabrics and dose it according to the manufacturer's instructions.
6. Please note that excessive washing will wear out the heating pad. For this reason, the heating pad should only be washed a maximum of 10 times throughout its entire life.
7. Only use the appliance once it is complete dry.
8. Do not switch the appliance on for drying!

### **Additional information:**

- The appliance can be used as a heating pillow or used as a heating pad after taking out the inside cushion.
- The SISSEL® heating pillow has been specially designed as warming cushion, not as sleep pillow, do not use it as sleep pillow.

### **Technical data:**

Size: approx. 35x46 cm Voltage: 230V~ 50Hz Watts: 100W
--

### **Disposal:**

The SISSEL® heating pillow must be disposed of in accordance with domestically applicable provisions.

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**SISSEL® GmbH**